

Zahnbehandlungsstand

Dental Treatment Table incorporating Mouth Gag



Zahnbehandlungsstand für
Nagetiere und Kaninchen das
Original

Mouth Gag
for rodents and rabbits
the original

Zahnbehandlungsstand (Maulöffner)

Der Zahnbehandlungsstand ermöglicht ein sicheres Öffnen des Mauls bei nahezu allen Arten von Nagetieren und Kaninchen. Die umfassenden Verstellmöglichkeiten gewährleisten die Anpassung an die jeweiligen anatomischen Verhältnisse und ermöglichen ein bequemes Arbeiten. Bei sedierten Tieren ist die Zahnkorrektur problemlos auch von einer Einzelperson durchführbar, da das Tier in der optimalen Position gelagert und gleichzeitig fixiert wird.

Der gleichmäßige Druck auf Ober- und Unterkiefer minimiert das Risiko von Verletzungen. Auch bei kurzen, fehlenden oder fehlgestellten Nagezähnen ist ein problemlos Öffnen des Mauls möglich, da der Mundöffner hinter den Nagezähnen ansetzt. Gleichzeitig wird durch diese Technik ein besserer Überblick in der Maulhöhle erzielt. Auch das Einsetzen des Wangenspreizers und die Arbeit mit den anderen Zahnbehandlungsinstrumenten werden deutlich erleichtert.


Rodent mouth gag (Mouth opener)

The rodent mouth gag provides a safe means of holding open the jaws of virtually all rodent and rabbit species. The wide range of adjustments ensures that the instrument can be adapted to the particular anatomical conditions, while also enabling the veterinary surgeon to adopt a comfortable posture. Tooth trimming is easy to perform on sedated animals even for a single operator because the animal is secured in the optimal position.

Simultaneous pressure on the upper and lower jaw minimises the risk of injury. Even in rodents with short, missing or malpositioned teeth, the jaws are easily opened, since the mouth opener is positioned behind the teeth. At the same time, this technique gives a better view of the oral cavity and also makes it easier to use the cheek dilator and other instruments for the dental treatment of rodents.

Lagerungsbeispiele

Storage example

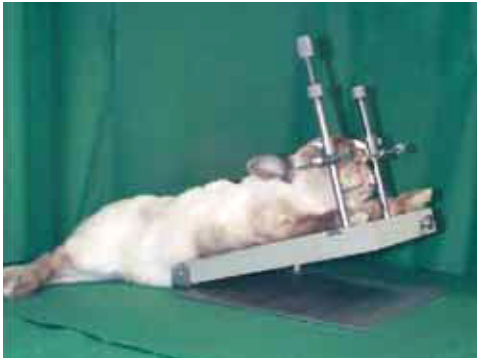

Lagerung eines Kaninchen
Positioning of a rabbit



Anwendung der Korrekturzange für Backenzähne
Application of adjustment forceps for molars



Abfeilen von Zahnspitzen
Using the dental rasp to file of long tips



Maulöffnen beim Chinchilla
Mouthgag in place (chinchilla)